



A family of high-performance crossing catheters for all of your peripheral interventions

Purpose-designed for peripheral treatment via radial access with lengths available for femoral and alternate access


The Sublime™ Microcatheter features our best-in-class Pristyne™ hydrophilic coating for superb lubricity. The proprietary shaft technology includes dual stainless steel, locked-braid layers which work in unison to offer enhanced torque control, push transmission and deliverability when tracking to distal target lesions. The PTFE-lined inner lumen provides smooth device tracking over guidewires and when telescoping microcatheters.




40 cm length Pristyne™ hydrophilic coating designed to optimize device trackability through vessel tortuosity

Working Lengths: 65, 90, 135, 150 and 200 cm

Guidewire Compatibility: .014", .018", .035"

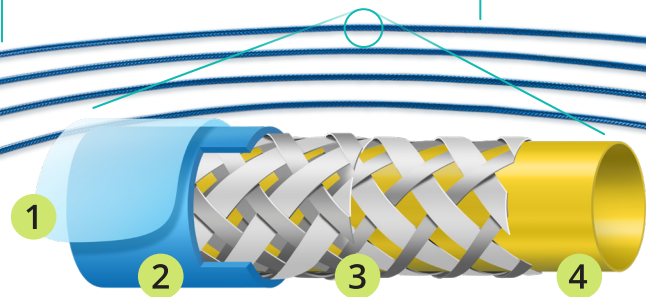

- 1 Pristyne™ hydrophilic coating designed to minimize vessel damage and spasm while optimizing trackability through distal tortuosity
- 2 Polymer outer jacket provides uniform device structure
- 3 Proprietary dual stainless steel, locked-braid layers work in unison to offer maximum torque control, push transmission, and deliverability
- 4 PTFE-lined inner lumen provides smooth device tracking over guidewires and when telescoping microcatheters

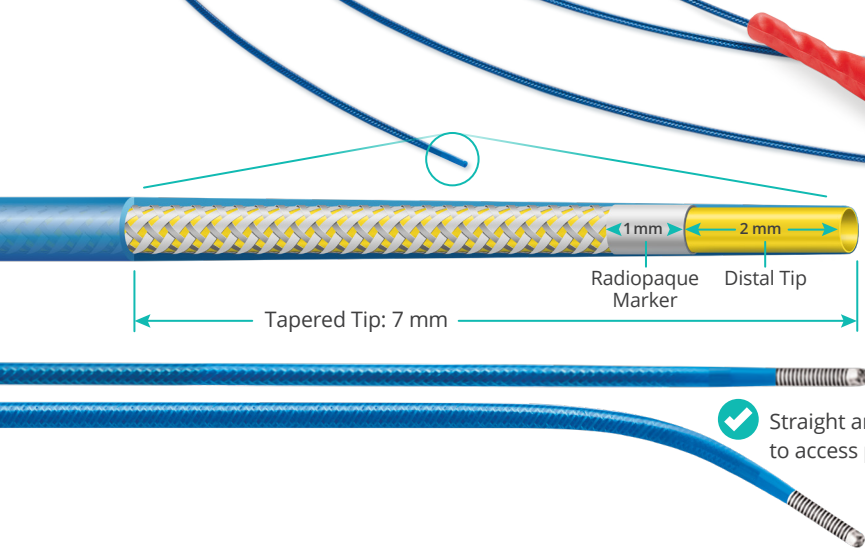


Straight and angled tip configurations to access peripheral vasculature




Sublime™ Microcatheter

The **.014"** and **.018"** microcatheter platforms can be telescoped through the **.035"** microcatheter platform to provide additional backup support when navigating heavily stenosed lesions.

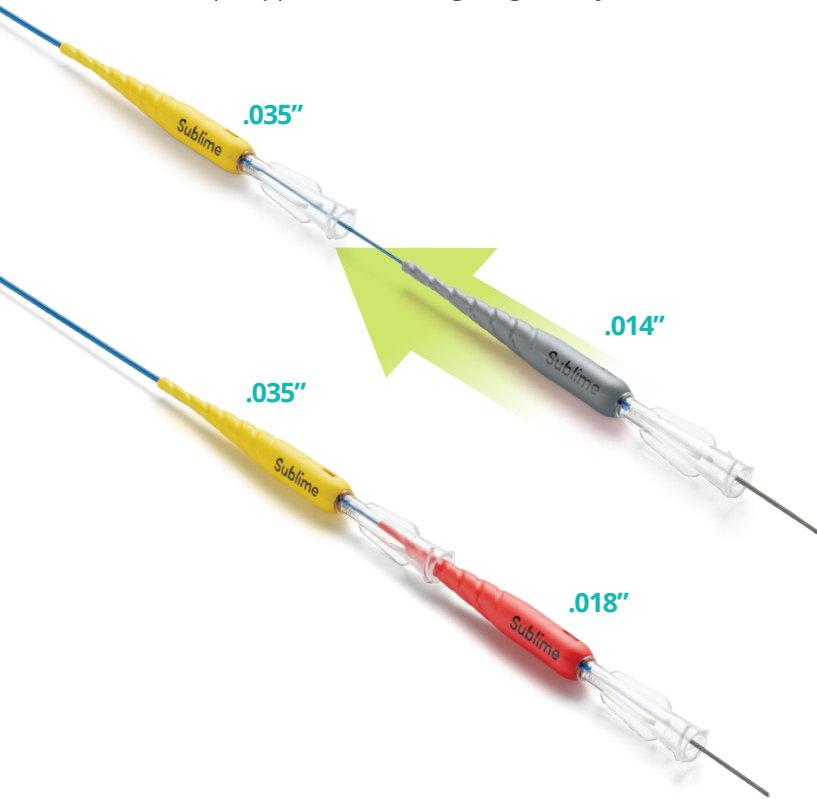

Model Number	Guidewire Compatibility (in)	Catheter Length (cm)	Tip Configuration
SRA-MC14-STR065	.014	65	Straight
SRA-MC14-ANG065	.014	65	Angled
SRA-MC14-STR150	.014	150	Straight
SRA-MC14-ANG150	.014	150	Angled
SRA-MC14-STR200	.014	200	Straight
SRA-MC14-ANG200	.014	200	Angled
SRA-MC18-STR065	.018	65	Straight
SRA-MC18-ANG065	.018	65	Angled
SRA-MC18-STR150	.018	150	Straight
SRA-MC18-ANG150	.018	150	Angled
SRA-MC18-STR200	.018	200	Straight
SRA-MC18-ANG200	.018	200	Angled
SRA-MC35-STR090	.035	90	Straight
SRA-MC35-ANG090	.035	90	Angled
SRA-MC35-STR135	.035	135	Straight
SRA-MC35-ANG135	.035	135	Angled
SRA-MC35-STR150	.035	150	Straight
SRA-MC35-ANG150	.035	150	Angled
SRA-MC35-STR200	.035	200	Straight
SRA-MC35-ANG200	.035	200	Angled

Working Lengths: **65, 90, 135, 150, and 200 cm**
 Guidewire Compatibility: **.014", .018", and .035"**

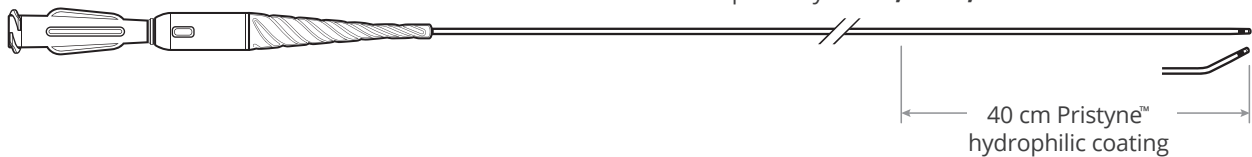

For additional information on the Sublime Microcatheter, please contact us at: contact@surmodics.com

INDICATIONS FOR USE

The Sublime Microcatheter is intended to access the peripheral vasculature in order to facilitate the placement and/or the exchange of guidewires. The Sublime Microcatheter is also intended to provide a conduit for the delivery of saline solutions or diagnostic contrast agents.

Caution: Federal (US) law restricts this device to sale by or on the order of a physician. Please refer to Instructions for Use for indications, contraindications, warnings, and precautions.



SURMODICS
 Surmodics, Inc.
 7905 Golden Triangle Drive
 Suite 190
 Eden Prairie, MN 55344 USA
 Toll Free: 888-626-8501
 Phone: 952-500-7400
www.sublimeradial.com

SURMODICS, PRISTYNE, SUBLIME and the SURMODICS and SUBLIME logos are trademarks of Surmodics Inc. and/or its affiliates.

©2022 Surmodics, Inc. All rights reserved. SRDX-SRA-710 A 11/22